



EZ9Plus & EZ11Plus User Self Help: “Door will not open”





Problem:

The power is turned off.
*(When the power is off
the door will lock)*

What to do:

Turn power on






Problem:

Sleep mode is active and the screen saver is displayed
(In the sleep mode the door will be locked)




What to do:

Awaken the unit from the sleep mode by pressing the START/STOP key.

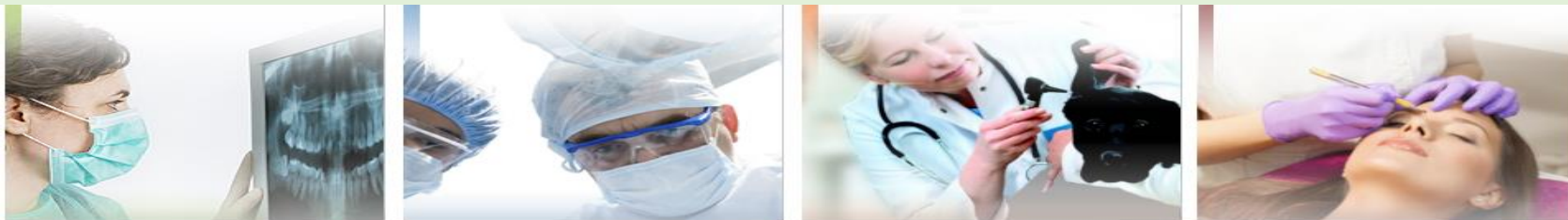



Problem:

There is water in the chamber and the water in the chamber icon  is displayed



What to do:

Start a new cycle. At the end of the new cycle water will be returned to the reservoir





Problem:
The machine has not returned to normal pressure

Wrapped Instruments, Pouches	
Ster. Temp.	270.0 °F
Ster. Time	4.0 min
Dry Time	30 min
System Ready	
Temperature	109.6 °F
Pressure	00.09 inHg

What to do:
Wait for the machine to return to normal pressure.



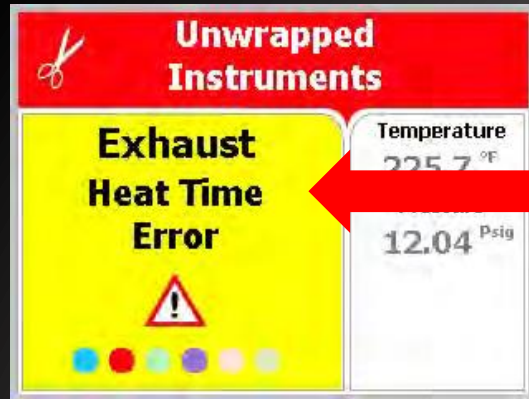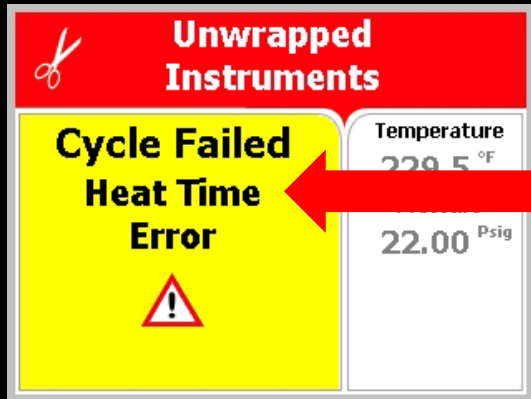

Problem:
An error message is displayed



What to do:
Press the START/STOP key to clear the error message and unlock the door

If the problem persists call your dealer or Tuttnauer USA for service.

